

OLYMPUS[®]

Your Vision, Our Future

ENDOALPHA

ENDOALPHA VIDEO MANAGEMENT

VIDEO NETWORK SOLUTION

The Smarter Way of Video Collaboration, Mentoring and Teaching




TOTAL VIDEO MANAGEMENT THROUGHOUT THE ENTIRE HOSPITAL

The ENDOALPHA Video Network Solution* (VNS) distributes and provides live videos wherever needed. Live or stored videos can be readily displayed everywhere within and beyond the hospital creating an ideal environment for teaching, mentoring and collaboration.

The total VNS package includes:

- **VNS web-based viewer (InStream Player™)** – completely web based player which can be accessed from any PC in the hospital.
- **Multi-Stream-Viewer** – enables multiple images to be viewed simultaneously on one monitor.
- **Session Manager** – coordinates the video network for the entire hospital.
- **Recording Suite** – allows for real time recording of any video source available on the network for teaching purposes and/or patient education as well as easy access to recorded content via the VNS web-based viewer.

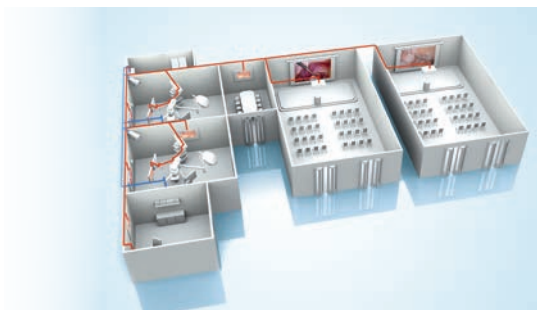

Features


- Allows hospitals to stream live video from the procedure room to any PC within and beyond the hospital network.
- Ability to view multiple streams (i.e. endoscopic, X-ray, pathology, room camera) simultaneously on a single PC monitor.
- Integrated access rights management, stand-alone or linked to LDAP or Active Directory.
- Full H.264 HD (1080p) and SD video support and full 256 bit video encryption.
- Instant and scheduled recording possible. Still image capture from video displayed on PC.
- Customisable screen layout for viewing preferences in Multi-Stream-Viewer.



Procedure / OR status:
Overview screen in the OR manager room

Benefits

- **No need** for doctors to **scrub-in** to mentor residents and students during procedures.
- Crystal clear HD images at doctor's office promotes convenient **collaboration, mentoring and teaching**.
- Remote participation for medical students during live procedures, therefore **reducing room traffic**.
- OR Manager is able to monitor **multiple room statuses** on a single screen, facilitating planning and tracking of patients and resources.
- Physician can easily reference and recall important **procedure events** by placing **HotMarks™** – video bookmarks embedded with meta data (recording).
- Can view procedure from any PC in the hospital network without installing specialised software or purchasing any additional licenses.
- **IT friendly** due to its web-based, installation free, zero footprint player.
- Remote service support allows **quick response times**.

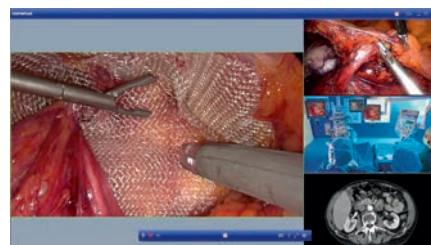

Screen view in the lecture room or doctor's office

* VNS may not be used for diagnostic purposes.

Olympus and ENDOALPHA are trademarks or registered trademarks of Olympus Corporation, Olympus America Inc., and/or their affiliates. Specifications, design, and accessories are subject to change without any notice or obligation on the part of the manufacturer.

OLYMPUS

OLYMPUS SINGAPORE PTE LTD. 491B, River Valley Road #12-01/04 Valley Point Office Tower, Singapore 248373, Singapore
OLYMPUS KOREA CO., LTD. Olympus-Tower, 114-9 Samseong-Dong, Gangnam-Gu, Seoul, 135-090 Korea
OLYMPUS (THAILAND) CO., LTD. 159 Serm-Mit Tower, 9th Floor, Unit #906-910, Sukhumvit 21 Road, North Klongtoey, Wattana, Bangkok 10110, Thailand

OLYMPUS (BEIJING) SALES & SERVICE, CO., LTD. A8F, Ping An International Financial Center, No. 1-3, Xinyuan South Road, Chaoyang District, Beijing, 100027 P.R.C.

OLYMPUS HONG KONG AND CHINA LIMITED L43, Office Tower, Langham Place, 8 Argyle Street, Mongkok, Kowloon, Hong Kong, China
OLYMPUS MEDICAL SYSTEMS INDIA PRIVATE LIMITED 102-B Time Tower M.G. Road Gurgaon 122001 Haryana., India
OLYMPUS AUSTRALIA PTY. LTD. 31 Gilby Road, Mt. Waverly VIC 3149, Australia
OLYMPUS NEW ZEALAND LIMITED 15d Paul Matthews Road, Albany, Auckland, New Zealand 1330